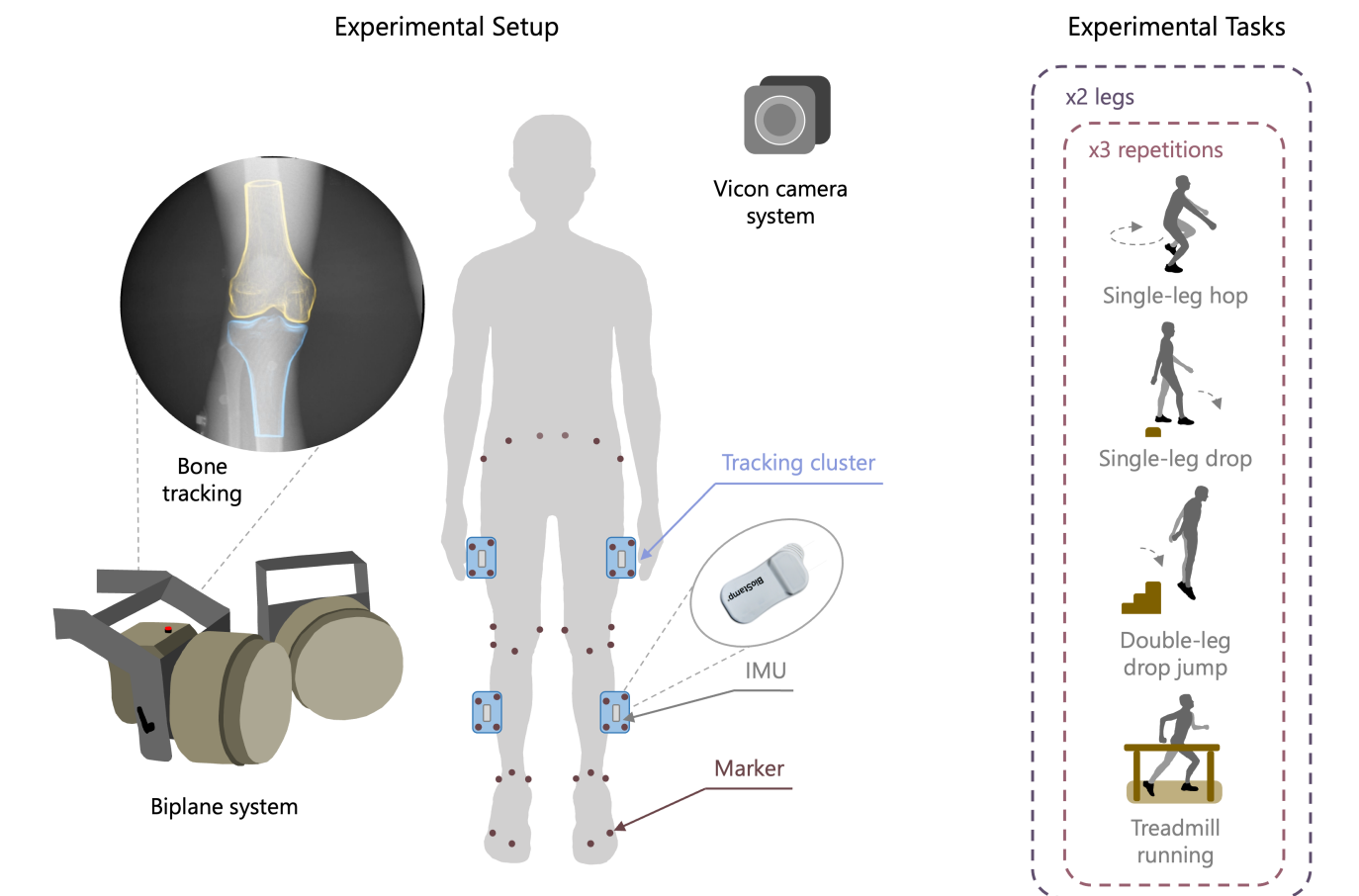




Experiment 4 (with Ground-Truth Biplane Fluoroscopy)

This folder contains the implementation to process data of Experiment 4.

Data

After downloading data from FigShare, follow the directory tree in [imu_biodynamics/data/README.md](#) to run the code.

Run the Code

Run `main_mc10.py` to obtain IMU kinematics.

```
python main_mc10.py
```

Run `main_mocap.py` to obtain marker-based kinematics.

```
python main_mocap.py
```

Run `eval.py` to obtain evaluation of IMU and marker-based kinematics against ground-truth biplane fluoroscopy.

```
python eval.py
```

Outputs from running these scripts are saved in the `outputs` folder.

Visualization

See `README.md` in the `outputs` folder to ensure to ensure that the data are formatted appropriately before plotting them.

If you do not wish obtain `outputs` from running the scripts above, you can also download the folder from FigShare and format it following the mentioned `README.md` file for visualization.

Run `plot_benchmark_f7_biplane.py` to plot the figure.

```
python plot_benchmark_f7_biplane.py
```